



Pass-Through Refrigerator for Lab/Clean Room

590

MODEL TEMPERATURE CAPACITY +15°C ÷ +2°C PT 5010PT

RANGE (LITRES)

 $(W \times D \times H \text{ cm})$ 67 x 89 x 204

EXTERNAL DIMENSIONS

DOOR TYPE

Dual F/R Glass Doors



OPTION: ESEDL-PT-2D: Electronic Lock with opening button on door handle





DRY CONTACT – Back-up Battery36hr - ACCURACY 0.1°C GAS R290 AVAILABLE

Colour: Gray, White or any RAL colour Structure: Internal and external Steel Laminated Plate with FINGER PRINT PROOF TECHNOLOGY™ to assure a Bacteria-Static Activity Insulation: Polyurethane Foam, CFC FREE, 42Kg/m3, thick. 60mm Refrigeration type: Ventilated (NO-FROST), Fast Recovery **Defrosting:** Automatic with condensing water evaporation, electric. Electric power: 1 / 210-240Volt / 50-60Hz

Power Protection: Shortage Protection, Power Stabilizer with Linear Voltage Regulator for power fluctuations (+/-20 Volt) Compressor: No.1 Hermetic Type - 600Watt - 3/4 Hp Refrigerant gas: R452a (Compatible with R404a) Door: No2 Glass Door with Electric Heater and Silicon Magnetic Gasket Double Fold, Self-Closing System with stopper at 95° Inner setting: 5 to 8 Shelves/Grids/Drawers Lights: 2 LED Light 120cm, 6000K

Electronic controller: Digital controller LED Further Display 7 segments microprocessor controlled, Accuracy 0.1°C - Data Logging Systems are available as option.

Additional READ OUT Temperature Controller on rear side Castors and Wheels: 4 Swivel Wheels + 4 adjustable feet

Main Visual Acoustic Alarm for: - HI and LOW temperature

- Open door
- Devices failure
- Probes failures

MAIN OPTIONS and ACCESSORIES:

DTLSD4--1T - Interactive Data Logger TFT 7" with USB/SD, Predictive Diagnosis, CAN-BUS Communication and faults tracking system, software for data read-out and graphics. etc. FLXP - Digital Data Logger with dedicated APP and WI-FI Connection, direct WI-FI data download and monitoring RG1 – Graphic Chart Recorder with 52 weekly Paper Disk and 1 Spare Pen.

TC - Twin Compressor, alternate work cycles, auto switch over if primary compressor is faulty

ESEDL-PT-2D – Electronic Lock with Password (Access Log) HLO - Hole for External Probe Ø 25mm

CS2 + DSR - Sliding Aluminum or S/Steel Drawers on bearings, edge 10cm H, Max Load 55Kg, Double Sliding System SH2 + DSR - Sliding INOX S/S perforated/plane Shelves on bearings Max Load 55Kg, Double Sliding System


SMS - SMS Alarm System in case of errors. WI-FI/3G - Online connection for Alarm Messages and Real Time monitoring, Graphics and Report Download GR - Remote Cooling System (max 8mts far from the unit,

version WIDER)

TR - Tropical Cooling System, for temp. till 45°C and 95% ru **NFOG** – NO FOGG Glass Technology

RS485 – Serial Port for PC/BMS Connection (only with DTLSD3 or FLXP)

TECHNICAL DRAWING with Standard Set-up

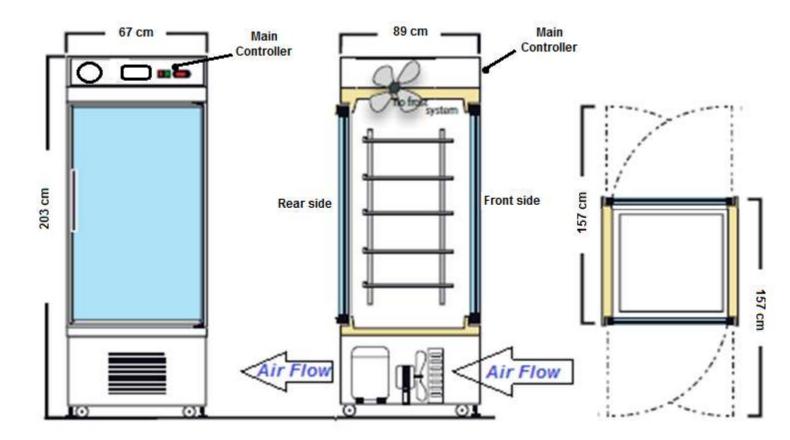

Options and Accessories

Sliding Drawers with Rear and Front pull-out rails (CS1-CS2 + DSR)



Temperature controller for static and ventilated refrigeration units at positive and negative temperature dedicated to the health sector; built-in WiFi module and Data Logger.



EMPLOYMENT UNIT

Range [-86°C ... +15°C]

- BLOOD BANK
- HOSPITAL REFRIGERATORS
- PHARMACEUTICAL REFRIGERATORS
- LABORATORY REFRIGERATORS

FUNCTIONS

- Monitor probe and the simultaneous display
- Test circuit and battery recharge
- Warning and recording of network failure
- Electric-key output driver
- Alarms and faults list with circular register
- Data Logger with clock and maintenance supercap
- Anti-ice function with safety thermostat
- Consistency check between the installed probes
- Buzzer with auto-silence and reset function
- Wireless parameters programming
- OEM Cloud supervision
- User supervision through app and web pages

1 Net

CE

- Weekly graphic report by mail
- Email alert messaging

PANEL DISPLAY

- Dimensions: L182 x H38 x D23mm
- Drilling template: L150 x H31mm
- Low voltage powerd by controller board
- Recessed, flush, under glass mounting
- 5 keys capacitive keyboard
- White and blue 7-segments display
- Data Logger (3Mb flash memory)
- Real Time Clock with super capacitor
- WiFi module built-in and alarm buzzer

CONTROLLER BOARD

- Power Supply: 230Vac 50/60Hz
- Nr 3 NTC probes [-49.9, +99.9]°C Res. 0.1°C
- Nr 2 PT1000 probes [99.9, +99.9]°C Res. 0.1°C
- Nr 4 Digital inputs (door switch, HP, LP, Thermal)
- Nr 1 Backup battery 12Vdc input
- Nr 6 Relais: 30A-10A SPDT-10A-10A-8A-3A
- Nr 1 Alarm relay SPDT
- Nr 1 Led bar 12/24V power supply input
- Nr 1 Led bar dimmered output [max 2A]
- Nr 1 Electric key output driver
- On-board connector for inserting modules: Ethernet, Modbus TCP, Modbus RS485, Inverter with electronic fan modulation outputs.





ap. Soc. 98.000 € 1 x - REA n. 218495 C.C.I.A.A. Varese. Reg. Impr. Varese. VA026-26187 - Cod. Fisc. e Part. IVA IT 01936060

Six excellent reasons to choose TEKNALAB

2

1

INTERACTIVE TOUCH PAD

A touch pad with 5 dynamic backlit keys simplifies navigation for the user. The central key is backlit by a RGB led which varies in colour according to the refrigeration unit operation.


Even from a distance it is easy to check the operating status of the refrigerator.

and to reduce consumption compared

Thermostat control with evaporator

probe in case of thermostat probe fail.

to traditional defrosting by time.

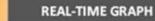

The temperature is shown on the screen with a graphical display that substitutes the now obsolete electro-mechanical thermograph.

The feature is integrated with a 10 year memory of the recorded data and list of events regarding temperatures and operation variables.


3

PREDICTIVE DIAGNOSTICS

In the mid-long term time period a relay fault is one of the most common failure to occur. In order to prevent this all relay switching are read with a representation of the remaining life expectancy of the component, and a warning is given when the threshold is about to be reached.



KEYBORD TOUCH WITH RGB BACKLIT

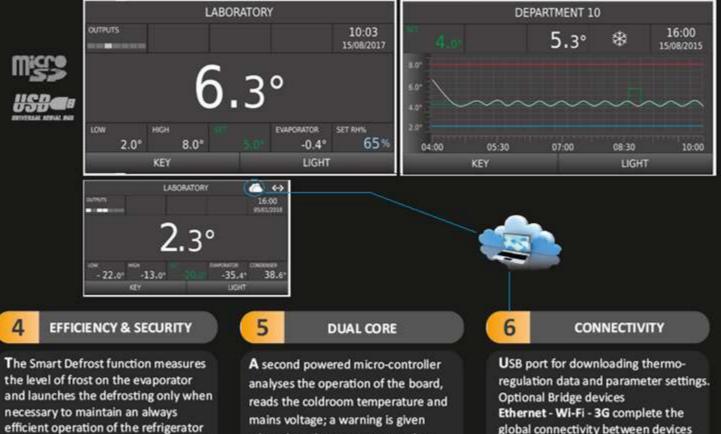

when the value are over or under

to safeguard compressor life and

electronics components.

the threshold defined by parameters

Ethernet - Wi-Fi - 3G complete the global connectivity between devices towards CLOUD Velex for monitoring and teleservice features via browser (PC or Smartphone) through credential access.

Control functions dedicated to each specific field of refrigeration: blood banks, laboratory and pharmacies

FUNCTIONS

Thermostat control

Regulation on central set with double hysteresis in cooling action for systems with single and double stage (superfreezers) Heating action for climatised rooms.

Defrosting

Defrosting mode: electric, hot gas, stop compressor with activation from the touch pad, by a timer, by a clock and by the automatic detection of ice.

Fan management


Condenser fans control within a range of temperatures to guarantee the correct thermal exchange.

Evaporator fans control to maintain humidity and to block the flow of hot air into the coldroom during defrost phase.

Programmable outputs

Any of the following actions can be associated to every relay: compressor, defrosting, evaporator fans, condenser fans, door resistance, glass resistance, drainage resistance, active load light, heating action and 2nd stage compressor action for superfreezers.

outputs wiring example

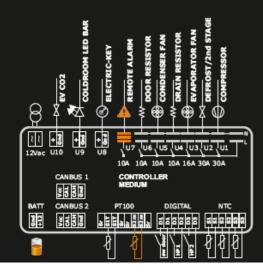

Redundancy

Double power supply for the logic.

A second micro-controller is designated to the function of diagnosis, battery charge and the reading and registration of temperature in normal conditions and in mains failure status.

For ventilated units in case of thermostat probe failure the regulation is transferred to the evaporator probe.

The temperature of the technical compartment is constantly monitored with a warning when it reaches the safety threshold.

The output and operation status is shown by exterior led (*expert points*).

PT100 probe monitor

PT100 precise independent monitor probe with a 0.1°C of resolution for thermo-registration and high and low temperature alarms with a coherence check between its readings and thermostat probe with alert for any possible unbalance.

DEDICATED FUNCTIONS

NTC/PT100 thermostat probe

Two thermostat probe inputs are provided to regulate the temperature both of the ventilated unit with NTC probe and of the superfreezers -86°C with PT100 probe.

Back-up battery

Back-up batteries are connected in a recharging circuit, with periodic tests, charging status and substitution alert. In case of a mains failure they provide 24 hours of power supply.

Electric-Key Driver

Solid state output controls an electric-key with password access.

Led bar driver

A solid state output can be configured to directly control the coldroom light by led bars with automatic functions of switching on, turning off and energy savings.

CO2 valve driver

A solid state output controls the CO2 solenoid valve in normal conditions and in the event of mains failure for the freezer units that foresee this.

Driver Board Power Supply 12Vac 50/60Hz >=6VA Dimensions L200 x H110 x P35 mm Inputs N°3 NTC [-50.0, 70.0]° N°3 Digital N°2 PT100 [-99.9, 50.0]° **Relais Out** N°2 x 30A + N°1 x 16A N°3 x 10A + Alarm 10A Drivers Led Bar - EV CO2 Electric-key Gates 2xCANBUS

Display Board

Power Supply	from Driver Board
Housing	plastic ABS
Dimensions	L160 x W52 x H25 mm
Dima	L154 x W49 mm
Keyboard	5 Keys Touch
	with RGB backlit
Display	TFT 2.41" 320x240 pixel
Gateway	CANBUS
	USB
	for data exchange