



Crucial Temperature Solutions

Twin *Guard*

-86°C Ultra-Low Temperature Freezers





575 L / 715 L

PHCbi's TwinGuard -86°C chest freezers deliver the unsurpassed level of safety and sample security that is critical for valuable biological sample storage alongside exceptional ease of use and data monitoring.

Ultimate Sample Protection

MDF-DC700VX

MDF-DC700VXC

The Dual Cooling System "TwinGuard" comprised of independent dual cooling circuits delivers -86°C, and one side of the cooling circuits maintains a chamber temperature at -65°C. The freezer minimises the risk of compromising valuable samples.

Improved Temperature Performance

TwinGuard chest freezers maintain better temperature stability and uniformity during routine door openings than equivalent upright models as less cold air escapes the freezer while the door is open.


Enhanced Use & Intelligent Security

A colour LCD touch panel allows full user control, even with gloved hands, while the USB port makes transferring logged data to a PC convenient. The filterless design of the freezers reduces routine maintenance.



Valuable Sample Storage

With two independent systems, each sustaining ultra-low temperatures, your samples remain safe. Optional CO_2 backup provides an extra layer of protection.

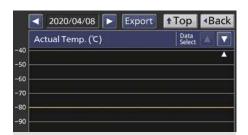

Predictive Maintenance

Performance monitoring ensures the freezer is operating within specifications. This protects the compressors and condenser fans. External stress on the compressors or fans, due to ambient temperature abnormality or excess load, is displayed for immediate notification.



Restricted Access Laboratories

Significantly extended critical time to react to any sudden mechanical failure. Ideal for high security applications.

TwinGuard

-86°C Ultra-Low Temperature Freezers



Filterless Design


The filterless construction of the freezers reduces routine maintenance time by eliminating the need for regular cleaning of filters.

Improved Temperature Performance


The TwinGuard -86°C chest freezer design provides optimal temperature stability and uniformity during routine lid openings, minimizing the loss of cold air even when the freezer lid is open.

Performance Data


AT30°C Pull-down & Pull-up Temperature

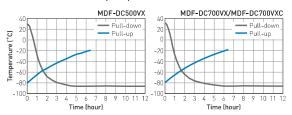

Dimensions

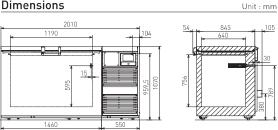

MDF-DC500VX

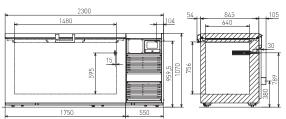

MDF-DC700VX MDF-DC700VXC

Model Number		MDF-DC500VX-PE	MDF-DC700VX-PE	MDF-DC700VXC-PA
External dimensions (W x D x H) 11	mm	2010 x 845 x 1070	2300 x 845 x 1070	
Internal dimensions (W x D x H)	mm	1190 x 640 x 756 1480 x 640 x 756		
Volume	litres	575 715		
Net weight	kg	328	358	366
Capacity	2" boxes	416	5:	20
Performance				
Cooling performance 21	°C	-86		
Temperature setting range	°C	-50 to -90		
Temperature control range 2)	°C	-50 to -86		
Control				
Controller		Microprocessor, non-volatile memory		
Display		LCD Touch Screen		
Temperature sensor		Pt-1000		
Refrigeration				
Refrigeration system		Independent Dual-Cooling		
Compressors	W	2 x 1100		
Refrigerant		HFC mixed		
Insulation material		PUF / VIP PLUS		
Construction				
Exterior material		Painted Steel		
Interior material		Stainless Steel		
Outer doors/lid	qty	1		
Outer door lock	- ' '	Y		
Inner doors/lid		3 pieces (Styrofoam)		
Access port	qty	1		
Access port position	.,	Back		
Access port diameter	Ømm	17		
Alarms (V = Visual Alarm, B = Buzzer Alarm, R = Remote Alarm)				
Power failure		V-B-R		
High temperature		V-B-R		
Low temperature		V-B-R		
Filter		Filterless design		
Door open		V-B		
Electrical and Noise Level		MDF-DC500VX-PE	MDF-DC700VX-PE	MDF-DC700VXC-PA
Power supply	V	230	/ 240	220
Frequency	Hz	5	0	60
Noise level 3)	dB (A)	52		
Options				
Liquid CO ₂ back-up		MDF-UB5-PW		
		MTR-G85C-PE		
စ္ - Circular type		— Chart paper: RP-G85-PW		
- Circular type		— Ink pen: PG-R-PW MTR-85H-PW		
S - Continuous strip type		мтк-85н-РW — Chart paper: RP-85-PW		
= 0 - continuous stub type		— Ink pen: DF-38FP-PW		
— Recorder housing: MDF-S3085-PW				
Optional Communication Systems			MTD LOG DIA:	
Ethernet interface (LAN) 4		MTR-L03-PW		
Digital interface (RS232C/RS485) 4)		MTR-480-PW		
Quality Management System				

- 1) External dimensions of main cabinet only, excluding handle and other external projections - See dimensions drawings on website for full details.
- ^{2]} Air temperature measured at freezer centre, ambient temperature +30°C, no load.
- 3) Nominal value Background noise 20dB[A]

Certification

^{4]} Only for MTR-5000 (data acquisition system) users.

Caution: PHC Corporation guarantees this product under certain warranty conditions. However, please note that PHC Corporation shall not be responsible for any loss or damage to the contents stored in the product.

IS09001



Preservation Equipment, Experimental Environment Equipment, Dispensary Equipment, Culturing Equipment and Drying & Sterilising Equipment for General Laboratory use

The management of the design, development, production and servicing of the above. PHC Corporation, Biomedical Division

1-1-1 Sakada, Oizumi-machi, Ora-gun, Gunma 370-0596, Japan





PHC Corporation Biomedical Division is certified for:

Environmental management system: IS014001

DISTRIBUTED BY:

